


Cari Chair

carilex®

Pump - Battery Powered/Portable

Giving you, the care you need whenever, wherever!

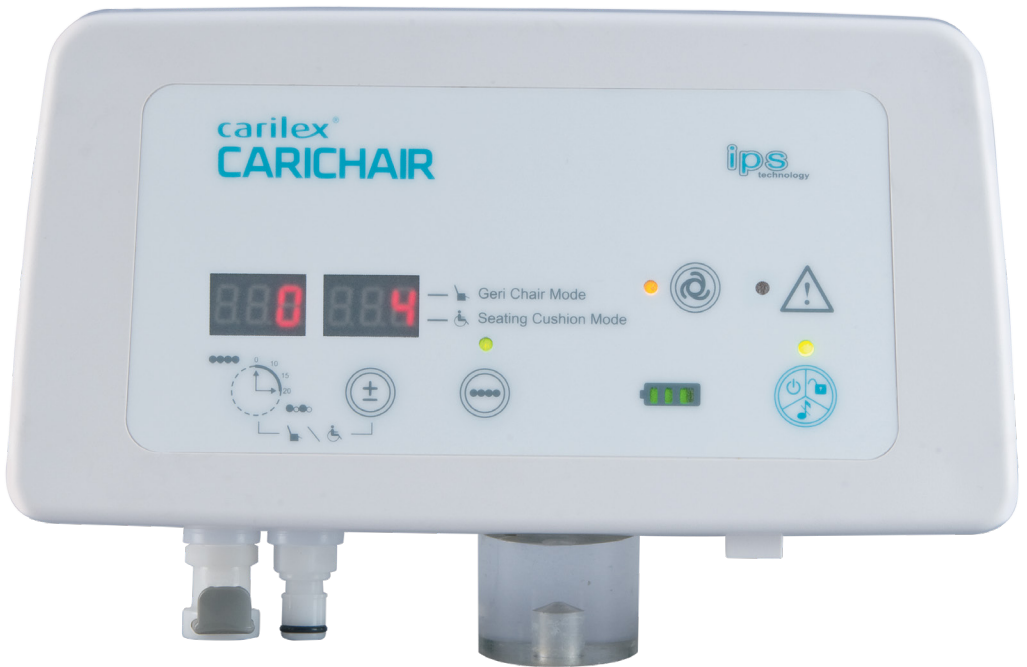

Whisper Quiet



Digital LED screen



Panel Lock



Visual and Audible indicator



Battery level LED



Alarm Mute

Global Leader in Wound Care, Distributed in Over 40 Countries Worldwide

Semi-automatic mobile cushion system with battery backup – defying the constraint of time and place for uninterrupted around the clock care. Whenever, wherever. Patients no longer need to be confined to the bed for pressure ulcer treatments, compatible with a variety of mobile devices from chair, recliner or wheelchair.

GENERAL FEATURES

Multiple Compatibility	Compatible with multiple cushions for various types of chairs
Intelligent Pressure Sensing	Responds to patient movements on cushions by automatically adjusting internal mattress pressure
Comfort Setting	Provides optimal internal pressure according to each patient's sensation. Includes options for comfort level adjustment, no weight input needed
Battery Backup up to 8 hours	Automatic power detection sensor activates battery backup, in the event of main power failure
Dynamic Mode	2:1 alternating therapy
Static Therapy/CLP	All air cells maintain constant pressure support
Multiple cycle times	10,15 and 20 minute cycles for alternating pressure therapy

TECHNICAL DETAILS

Risk/Category	Med - High
Waterlow Score	15+
Pressure Sore Grading	3
Safe Working Load	Up to 150kgs
Compatible Support Surfaces	Cari Chair, MobiCare and MobiCare Plus
Travel/Weight	Weighing 0.92kg
Clinical Settings	Suitable for Hospital and Home Care environments
Warranty	1 year

carilex®

Manufacturing for over 30 Years - 10 Years in the Australian Market