

Coziny Plus ^{carilex}®

Paediatric Pressure Care

Carilex Coziny Plus offers wound care therapy designed specifically to take care of children who are at risk of developing pressure ulcers, or already enduring distress and discomfort from pressure ulcers. Intelligent Pressure Sensing (IPS) Technology, responds to patient movements on mattress by automatically adjusting internal mattress pressure, which allows Coziny Plus to regulate airflow and continuously provide total envelopment for our patients.



Whisper Quiet



Memory Recall



Panel Lock



Alarm Mute




Battery level LED



Visual and Audible Indicator

Global Leader in Wound Care, Distributed in Over 40 Countries Worldwide

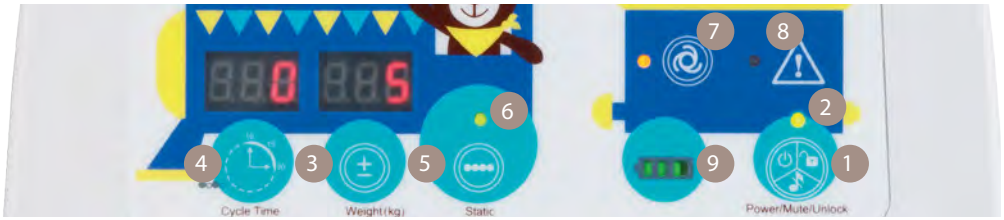

GENERAL FEATURES

1 Power/Mute/Panel Lock	Green LED indicates pump is on/Mute audible alarm/Hold for 3 sec to unlock pump
2 Power Failure Alarm	LED when red power failure. Audible and visual power failure alarms to alert caregivers in case of abnormal power outage
3 Multiple Weight Settings	Selectable weight settings to ensure optimal pressure relief for all body types
4 Multiple Cycle Times/Dynamic Mode	Various cycle time selection for alternating pressure therapy, 10, 15 and 20 minutes. Dynamic 2:1 alternating cycle
5 Static Therapy	Constant low pressure mode, all cells in the mattress
6 LED Light	Lights up green when active
7 Intelligent Pressure Sensing	Responds to patient movements on mattress by automatically adjusting internal mattress pressure
8 Low Pressure	Exclamation LED flashes Orange when insufficient internal pressure
9 Battery Backup up to 8 hours	Automatic power detection sensor activates battery backup in the event of main power failure and allows treatment to continue uninterrupted during patient transportation and examination
Transport Mode	Connect the two quick couplings together to retain air in the mattress, can stay inflated for up to 15 hours
Multiple Compatibility	Compatible with multiple mattress sizes for various body types

TECHNICAL DETAILS

Risk/Category	Medium to High
Waterlow Score	15+
Pressure Sore Grading	3
Safe Working Load	Up to 25kg
Compatible Support Surfaces	Coziny 100 (incubator), 200 (baby crib) and 300 (baby cot)
Travel/Weight	Pump weighs 0.92kg
Clinical Settings	Suitable for Hospital and Home Care environments
Warranty	2 year

carilex®

Manufacturing for over 30 Years - 10 Years in the Australian Market